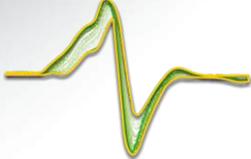


Version 7

Signal 

Sweep-based data acquisition & analysis system

Signal is a sweep-based data acquisition and analysis package. Uses range from a simple storage oscilloscope to complex applications requiring stimulus generation, data capture, control of external equipment and custom analysis. This flexibility makes it ideal for a wide range of applications including transient capture, patch and voltage clamp, LTP studies and evoked response.

- Signal is simple to set up for data capture and analysis
- Signal provides, in a user-friendly environment, the tools for sampling and analysis that most researchers need. A built-in script language automates repetitive tasks with additional functions for custom analyses and applications
- Signal includes functions for specific application areas including dynamic clamp, whole cell and patch clamp electrophysiology and evoked response studies with control of magnetic and other stimulus devices
- Signal imports data recorded by many other systems, so you can take advantage of this extremely versatile package to analyse existing data

The powerful data capture and analysis functions make Signal, together with one of the CED1401 family of laboratory interfaces, an extremely flexible, cost-effective and time-saving addition to any laboratory.

Typical applications

Evoked response and TMS Output fixed, random and pseudo-random sets of pulses with on- and off-line waveform averaging and measurements of latencies, amplitudes and areas. Signal can control supported transcranial magnetic stimulators during data acquisition, including adjusting stimulator amplitude and timing with checks on stimulator condition. Settings are stored in the corresponding data frame.

LTP and LTD studies Generate single, paired and trains of pulses on multiple outputs. Automatically measure features including baselines, amplitudes, latencies, areas, durations, percentage rise and decay times, and population spike parameters for single or multiple responses per frame.



Whole cell voltage and current clamp Interactively control holding potentials and pulse outputs using a graphical pulse editor. Construct I/V plots, leak subtraction and curve fits both on- and off-line.

Dynamic clamp High speed, fully integrated functions to simulate ion channels and leaks and use simulated synapses to generate hybrid networks.


Single channel patch clamp Detect single and multi-level channel openings on-line and off-line automatically. Produce measurements such as open/closed times and amplitudes and display result histograms.

EMG Record single and multiple channels of EMG with software control of amplifiers such as the CED1902. Perform rectification and smoothing at the touch of a button or under script control. Produce power spectra that update as new data sweeps are captured.


Auditory evoked potentials Generate complex stimulus protocols and reports under script control. Features include: artefact rejection, digital filtering, generation of sub and grand averages and feature detection using active cursors.

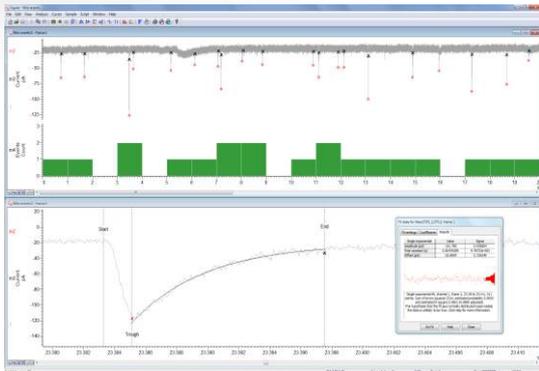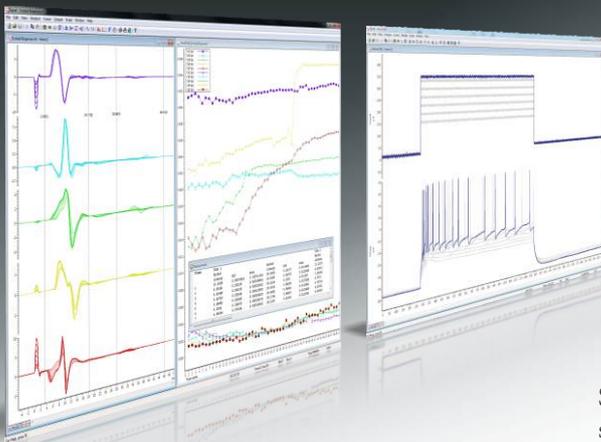

Automated detection and measurement of synaptic events

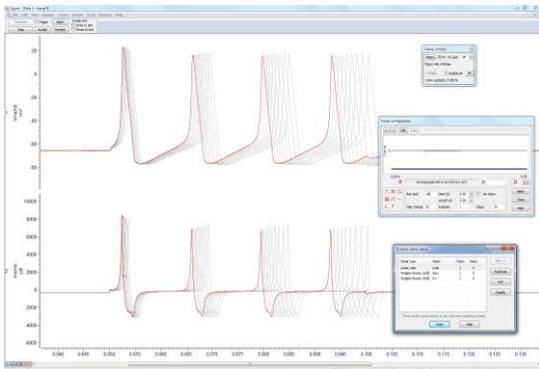

Action potentials simulated by dynamic clamp

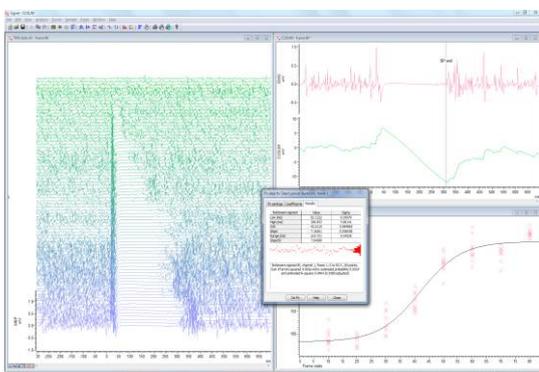

TMS: Measurement of silent period duration using CUSUM method

Data acquisition using a CED 1401


- Record and display fast sweeps of waveform data to disk with up to 40 million data points per sweep with a Power1401-3A or towards 10 million with a Micro1401-4
- Output pulse, sine, ramp and stored waveforms while sampling
- Generate sets of varying stimuli on up to 8 waveform and 16 digital outputs
- Set up multiple output states, each providing differing pulse outputs or stimuli. Sequence the order and repeats of each state cyclically, randomly, semi-randomly, manually or to a predefined protocol. Frames are marked with the state used
- Sample multiple waveform inputs at aggregate rates up to 1MHz
- Easily switch sampling configurations or output protocols
- Capture free-running or triggered sweeps with post or peri-trigger modes
- Gap free option for acquisition of continuous data
- Trigger data capture from external pulse, waveform threshold crossing or internal clock
- Sample variable numbers of points in different sweeps
- Tag frames to indicate state or artefacts for selective review and analysis
- Perform built-in and custom analyses on incoming sampled data
- Configure software controllable amplifiers including the CED1902 low noise isolated pre-amplifier, the Digitimer D360 and D440 isolated patient amplifiers and the Power1401 gain option
- Automatically adjust channel scaling in response to gain changes detected from either patch clamp amplifier telegraphs or supported software controlled amplifiers
- Simultaneously review earlier frames together with incoming data using duplicate windows
- Sample at arbitrary or varying rates with external convert input option

Complex stimuli and experiment control


Signal incorporates a powerful sequencer function for generating digital and analogue outputs for stimulus control. Sequences are run within the 1401 interface to ensure highly accurate timing independent of the host PC operating system. There are two methods available for specifying outputs.

The graphical pulse editor covers the majority of stimulus requirements in an easy to use drag and drop environment where you build up to 256 sets of up to 500 pulses linked to the multiple states system. Outputs are fixed or they can change amplitude and duration on repeats. You can also modify the outputs while sampling continues.

If your requirements cannot be entirely met with the pulse editor you can define your outputs as a text sequence of instructions. This makes it possible to control the sequence interactively or by using the script language running on the host PC to pass variables and tables of values to the 1401. It also enables you to produce pulse protocols either time-locked to the recording frames or free running. Further experiment control options including communication through serial lines are also available using the script language.

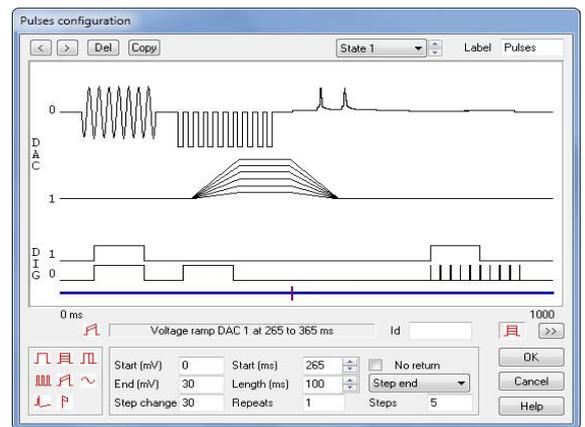


Output pulse editor

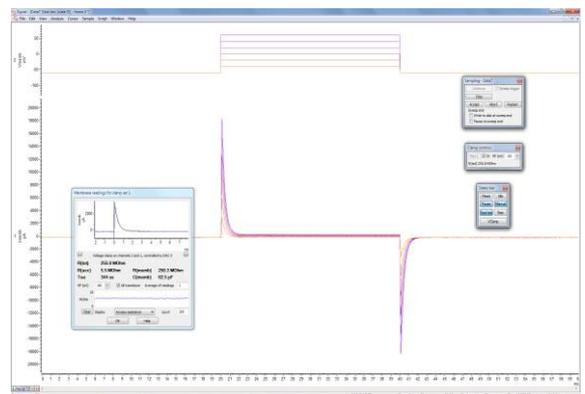

Voltage clamp recording using model cell



MEP recording system

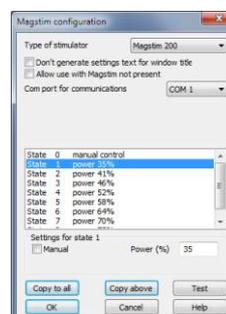

Magstim configuration dialog

State	Percentage RMT	Magstim level
1	60	35
2	70	41
3	80	46
4	90	52
5	100	58
6	110	64
7	120	70
8	130	75
9	1	1
10	1	1

Script control of Magstim power settings. This example sets power levels relative to RMT




CED 1902 Isolated pre-amplifier

Script language

From complete beginners to experienced programmers, everybody can benefit from the built-in script language. Even simple automation of repetitive tasks with known parameters can save hours or even days of tedious analysis. Advanced scripts can provide complete experiment control with on-line application of original algorithms to sampled data in real-time.

The script language allows you to generate your own user interface, and there is a macro recording facility to provide a starting point for new scripts. The script language also includes data manipulation tools such as multi-dimensional arrays and matrix functions.

CED maintains a library of example scripts for a wide range of common and specialised applications. If the scripts included with the Signal software or those available from our web site do not fulfil your exact requirements, please contact us to discuss the available options. These include modification of existing scripts as well as a dedicated script-writing service.

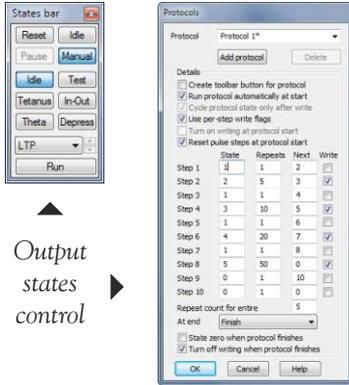


Output states control

System requirements

For sampling, Signal version 7 requires a CED Power1401, Micro1401 or 1401plus laboratory interface and a PC with Windows 7/8/10 or an Intel Macintosh running Windows. We recommend a minimum of 4GB of RAM. Signal dynamic clamping requires a CED Power1401 mk II or Power1401-3 laboratory interface.

Both 32-bit and 64-bit versions of Signal are supplied; either version of the software can be used on systems running a 64-bit OS. The 64-bit version provides a performance increase and access to more memory.

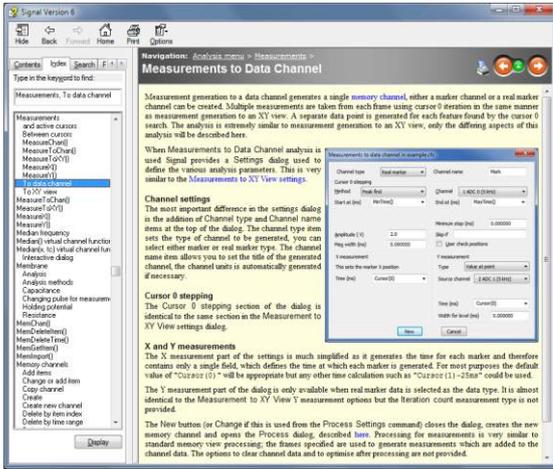

Signal on-line help



Signal user training days



www.ced.co.uk

Cambridge Electronic Design Limited

Technical Centre, 139 Cambridge Road, Milton, Cambridge CB24 6AZ, UK. **Tel:** (01223) 420186
Email: info@ced.co.uk **Europe & International Tel:** [44] (0)1223 420186 **USA and Canada Toll free:** 1-800-345-7794
Distributors in: Australia, Austria, China, France, Germany, Israel, Italy, Japan, Switzerland & Turkey